

## Product Information

<b>Name</b>	Anti tPAI-C monoclonal antibody
<b>Target</b>	tissue plasminogen activator–inhibitor complex, tPAI-C
<b>Target Species</b>	Human
<b>Target Description</b>	tPAI-C is the complex between t-PA and PAI-1, which reflects the severity of endothelial cell damage and the resulting level of fibrinolysis.
<b>Host</b>	Mouse
<b>Product Description</b>	Monoclonal antibody produced by <i>in vitro</i> cell culture under conditions free from animal-derived components
<b>Conjugation</b>	Unconjugated/AP/AE/HRP
<b>Catalog #</b>	1020201/1020201(AP/AE/HRP)
<b>Solution</b>	Saline solution, pH 6.5-7.5 with preservative
<b>Appearance</b>	Clear, transparent liquid, colorless or slightly yellow
<b>Storage</b>	2-8°C
<b>Purity</b>	>95%
<b>Subclass</b>	IgG1, kappa
<b>Affinity (K<sub>aff</sub>)</b>	>10 <sup>9</sup> L/mol
<b>LoB</b>	User Guide
<b>Specificity</b>	Positive: / Negative: TM, TAT, PIC
<b>Epitope</b>	ND
<b>Applications</b>	<input checked="" type="checkbox"/> ELISA <input checked="" type="checkbox"/> CLIA <input type="checkbox"/> GICA <input type="checkbox"/> FIA <input type="checkbox"/> IHC
<b>Pair</b>	<input type="checkbox"/> Capture <input checked="" type="checkbox"/> Detection
<b>Recommendations</b>	Pairable with #1020101
<b>Stability</b>	2 years
<b>Concentration</b>	2.0 mg/mL (+/-10%) or marked on the vial
<b>Lot #</b>	_____
<b>References</b>	

## User Guide

Catalog # 1020

Lot # \_\_\_\_\_

### CLIA

#### Antigen Identification


The antigen adopted in the following assays is a widely recognized clinical applied product obtained commercially.

#### Procedure

1. Mix 20 $\mu$ L serially diluted antigen or serum samples with 20 $\mu$ L (diluted to working conc. 50T/mL) capture reagent of 1020101(B) and 20 $\mu$ L 0.6mg/mL streptavidin-conjugated magnetic beads. Incubate with stir at 37 $^{\circ}$ C for 30 minutes.
2. Add 50 $\mu$ L (diluted to working conc. 20T/mL) detection reagent of 1020201(AP). Incubate with stir at 37 $^{\circ}$ C for 30 minutes.
3. Measure the relative light unit (RLU) of each reaction with a clinical approved chemiluminescence analyzer.

#### Measuring Range

The linear measuring range is 0.451-100ng/mL ( $R^2>0.99$ ).

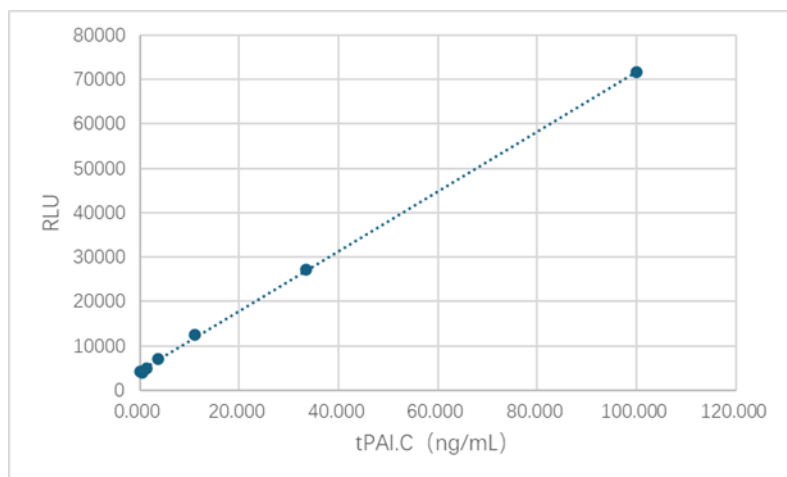

**Limit of Blank\*** \*Follow the Guideline EP17 published by Clinical and Laboratory Standards Institute (CLSI)

Sample tested: blank samples. Tests: 20.

The LoB is <0.5ng/mL.

**Limit of Detection\*** \*Follow the Guideline EP17 published by Clinical and Laboratory Standards Institute (CLSI)

Sample tested: low concentration samples. Tests: 20.

The LoD is 0.451ng/mL.

#### Specificity

The assay shows no cross-reactivities when detecting antigens as follows: TM at a tested conc. of 500ng/mL, TAT at a tested conc. of 500ng/mL, PIC at a tested conc. of 500ng/mL.

**HKIG™ All Rights Reserved**